medolight® BluDoc®

EN	Instruction manual	2
BG	Инструкция за употреба	8
CS	Použití přístroje	14
DE	Bedienungsanleitung	20
FR	Mode d'emploi	26
HR	Medolight upute za korištenje	32
HU	A készülék ki- és bekapcsolása	38
IT	Manuale d'uso	44
PL	Instrukcja obsługi	50
RO	Manual de utilizare	56
RU	Инструкции по использованию	62
SK	Návod na použitie	68
SR	Uputstvo za upotrebu	74
UA	Увімкнення та вимкнення апарата	80
AR	تعليمات الاستخدام	86

Security measures

- Do not make the device available to children – the device can only be used by adults who thoroughly read the manual.
- The AC adapter YS10MB should be connected only to a corrected installed mains socket. Make sure that the voltage in the mains socket is compatible with the voltage listed on the nameplate of the adapter.
- Arrange the cable connected to the device so as there is no possibility of accidentally tripping over it.
- In case of a liquid or foreign body getting into the housing of the unit, immediately turn the device off and send it to an authorised service centre for inspection.
- Do not use organic solvents for cleaning the housing. Clean the device with a soft, damp cloth. Before cleaning remove the AC adapter from the mains socket.
- In addition to the visible light, the device emits infrared (IR). While radiating with medolight devices avoid direct exposure of the eyes and close them when radiating areas close to the eyes.
- Wearers of contact lenses should remove them for the duration of the therapy and keep eyes closed.
- Do not immerse the AC adapter and the device in water and pour water over it, do not wash in the dishwasher. If necessary, wipe with a lightly moistened cloth
- Do not use the device showing any signs of damage: damage to the connection cable, housing and others.
- Device can only be repaired by authorised ZEPTEK service centre. The manufacturer shall exclude any liability in case of any unauthorised repairs to the device.
- Do not repair or modify the device by yourself – the damaged device should be returned to an authorised service centre.
- **This equipment may be used by children older than 7 years old and persons with reduced physical or mental capabilities, as well as persons with no experience and knowledge of the equipment if they are supervised or instructed on using the equipment in a safe manner so that the risks associated with its usage are understood. Children should not play with the product. Children without supervision should not carry out the cleaning and maintenance of the device..**
- Pregnant women should consult their doctor before the therapy.
- People suffering from acute or chronic diseases should consult their doctor whether this kind of physical treatment is recommended.
- PEOPLE SUFFERING FROM CANCER, PHOTOALLERGIES AND PREGNANT WOMEN SHOULD NOT USE THE medolight DEVICES. IN THE CASE OF EPILEPSY DO NOT USE THE DEVICE TO RADIATE HEAD NOR USE THE PULSE LIGHT. DO NOT USE THE DEVICE FOR RADIATING THYROID AFFECTED BY A DISEASE.
- The device is a source of heat

Description and characteristics of the device

The set consists of the medolight device and the AC adapter. The device is used to radiate with light emitted by LED in the range of visible light and infrared light wavelengths.

Principle of operation

The medolight BluDoc device is equipped with a matrix consisting of 171 blue LEDs (emitting light with wavelength of 470nm) and 36 infrared LEDs (emitting infrared light with a wavelength of 940nm). The device is powered up by a built-in battery which can be charged using the AC adapter included with the device. It is possible to use the device while the battery is being charged. However, it increases the charging process. The unit is controlled by a microprocessor. You can select one of five programs and one of five radiation durations. The processor controls the battery charging status and indicates the need of charging.

How to use

The medolight devices are an important part of a daily body care.

1. The device can be used 1-2 times a day for 1-2 cycles (1 cycle from 5 to 25 minutes). The device can be used at any time when you have time and opportunity to enjoy biostimulation. According to specialists the best results are achieved when using medolight device in the morning and evening.
2. Before using the medolight device, clean and disinfect the area subjected to radiation.
3. Relax and direct the stream of light to the place undergoing treatment from a distance of about 5 mm. While using the device, make sure that, if possible, the stream of light does not move. For large surfaces you should gradually direct the stream to next areas by moving the device to the non-radiated area after radiating the previous one.

Device Operation

1. Activation

Press and hold the ON/OFF button by at least 2s. Bottom icon bar (5-25) will indicate the battery status for 2 s – number of lit icons indicates the battery charging status. Then, program selection icon will light up (P1-P5) along with the time selection icon (5-25), in line with last settings before stand-by mode.

2. Program selection

Pressing the „P” button (PROGRAM) changes the radiation program signaled by lighting the appropriate icon with red light:

- P1 – continuous light
- P2 – pulsating light 10 Hz
- P3 – pulsating light 600 Hz
- P4 – pulsating light 3000 Hz
- P5 – pulsating light 8000 Hz

3. Setting the radiation time

By pressing the „T” button (TIME) you can select the radiation time in sequences 5 min → 10 min → 15 min → 20 min → 25 min → 5 min... signaled by the icon being backlit with green light.

4. Radiation

After selecting the radiating program and time, briefly press the ON/OFF button. LED matrix is now lit. Direct the light on the place to be radiated.


5. Pause

At any time, briefly pressing the ON/OFF button will turn off the LED matrix. Selection of a given program (P1-P5) and the time remaining to the end of radiation (5-25 min) are lit up.

If within 1 minute you briefly press the ON/OFF button again, radiating will be resumed. If within 1 minute you do not press the ON/OFF button again, the device will go into the sleep mode. During the pause it is possible to change the program and the radiation time by pressing the „P” (PROGRAM) or „T” (TIME) buttons (sections 2 and 3). You can end the radiation procedure by pressing the ON/OFF button for 2 seconds.

6. Radiation end

Radiating will automatically end after the selected time and the device will go into the sleep mode. Radiating can be finished by pressing the ON/OFF button for 2 seconds at any moment. After a short beep, the device will go into the sleep mode.

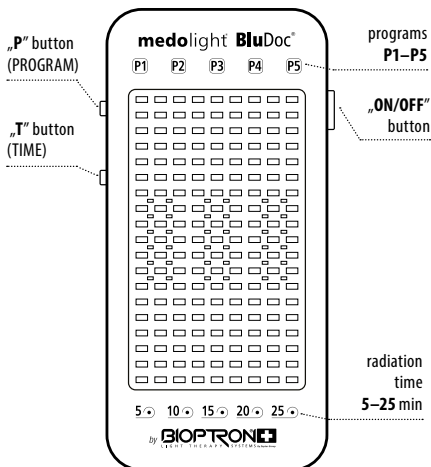7. Low battery indication

Low battery status is indicated by four short beeps. You must recharge the device – see section 8. If you are using the device, it will remind you that you need to recharge the device with four short beeps every 15s. If you do not plug the charger, the device will beep for 2 seconds and then go into the sleep mode.

8. Charging the battery

After you connect the cable and turn on the charger, some of the icons in the bottom bar will lit up proportionally to the battery voltage. Icons will be lighting up proportionally to the voltage increase. You will hear an acoustic signal for 1s when the charging is complete and the device will get deactivated. The last icon remains lit until you disconnect the charger. Turning the charger on during operation starts battery charging but is not indicated and does not interrupt the currently executed program function.

CAUTION: In the case of the battery being completely discharged, connecting the charger may not be immediately indicated. The indication may appear after a few minutes. It is not a fault but a natural feature of the battery.

9. Operating while recharging.

If you press and hold the ON/OFF button for 2s while the battery is recharging, the device gets activated. Normal device operation is possible. Battery recharging is independent from the radiating process. After the end of the radiation the device will start indicating the recharging until its completion.

Troubleshooting

Fault	Possible cause	Remedial measures
No response to pressing the ON/OFF button – the device cannot be activated	Discharged battery	Plug the AC adapter and charge the battery
	No voltage in the power outlet	Plug the AC adapter into a working socket
Device not responding when the AC adapter is plugged in	Faulty contact of the adapter plug	Unplug the AC adapter and plug it again by pressing well
	Faulty contact in the device socket	Unplug the USB plug from the device and plug it again by pressing well

Indications and symbols

present on the label and the Operating Manual

Operating Manual

Read the operator's manual

Symbol of the insulation protection class II

Manufacturer's declaration that the labelled product complies with the so called „New Approach“ directives of the European Union (EU) which apply to it

The EAC mark means that the product complies with applicable standards, has passed all Eurasian Customs Union conformity assessment procedures and is authorized for trading within the Customs Union countries.

Waste management – Proper disposal of this product. This designation indicates that in the EU the product should not be disposed with household waste. To prevent possible harmful influence of uncontrolled disposal of waste on the environment and human health, the waste should be disposed of in a manner allowing re-use of materials. In order to return your old appliance, please contact the services which collect such waste or the retailer where you purchased the product. Stores have the option to transfer the product to be recycled.

Warranty

All our devices have undergone very extensive quality control tests and we guarantee that no damage occurred to them in the production process. Warranty is valid for 2 years from the date of purchase indicated on the invoice issued by the seller at the time of delivery. The manufacturer undertakes to repair or replace elements in which structural defects were revealed during the warranty period. Faults which cannot be indisputably attributed to material defects will undergo expertise in our workshop and the cost of any repair will be priced depending on the result of the inspection.

The warranty is valid on condition that:

1. The device has material defects or faults not resulting from improper use.
2. The device does not have any traces of unauthorised repairs.
3. The device has not been damaged or broken as a result of user's negligence.

The warranty does not cover:

1. Completely destroyed parts.
2. Damage caused during transport, arising as a result of misuse, negligence or using the device in a way contrary to the warnings set out in the device manual.

The warranty does not affect the rights of the consumer, arising from the national legislation, nor their rights resulting from the agreement of sale.

The Warranty Card included with the product and the purchase invoice details terms of warranty.

Technical data

Power:

– external battery	3,7V 1450 mAh LI-ION
– YS10MB AC adapter	
– nominal feed voltage	100-240V, 50/60Hz
– nominal output voltage	5,4V ± 10%
– power consumption	max 12 W

Typical charging time 3 hours

Operating time after charging the battery for up to 60 minutes

The length of emitted light waves:

– infrared	940nm +/- 30nm
– blue	470nm +/- 30nm

The light energy in 5 min: max 2,75J/cm²

Exposure time up to 25 min +/- 5%

Ambient temperature limits:

– operation	+10°C ÷ +40°C
– transport and storage	- 20°C ÷ + 45°C

Dimensions of the device (L/W/H) 131 x 65 x12 mm

Weight:

– medolight device	112 g
– YS10MB AC adapter + cable	72 g
– case	11 g

Product complies with standards: EN 60335-1:2012 + A11:2014
 EN 60335-2-27:2013 + A1:2020
 +A2:2020
 EN 62233:2008 +AC:2008
 EN 62471:2008
 EN 55014-2:2015
 EN 61000-3-2:2014
 EN 61000-3-3:2013+A1:2019

Ver. 3 of 31.03.2022

The product meets the requirements of EU directives

Package contents:

- medolight device
- AC adapter
- Case
- Operating manual
- Warranty card

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